
BioRePeelCl$_3$


BIPHASIC TECHNOLOGY



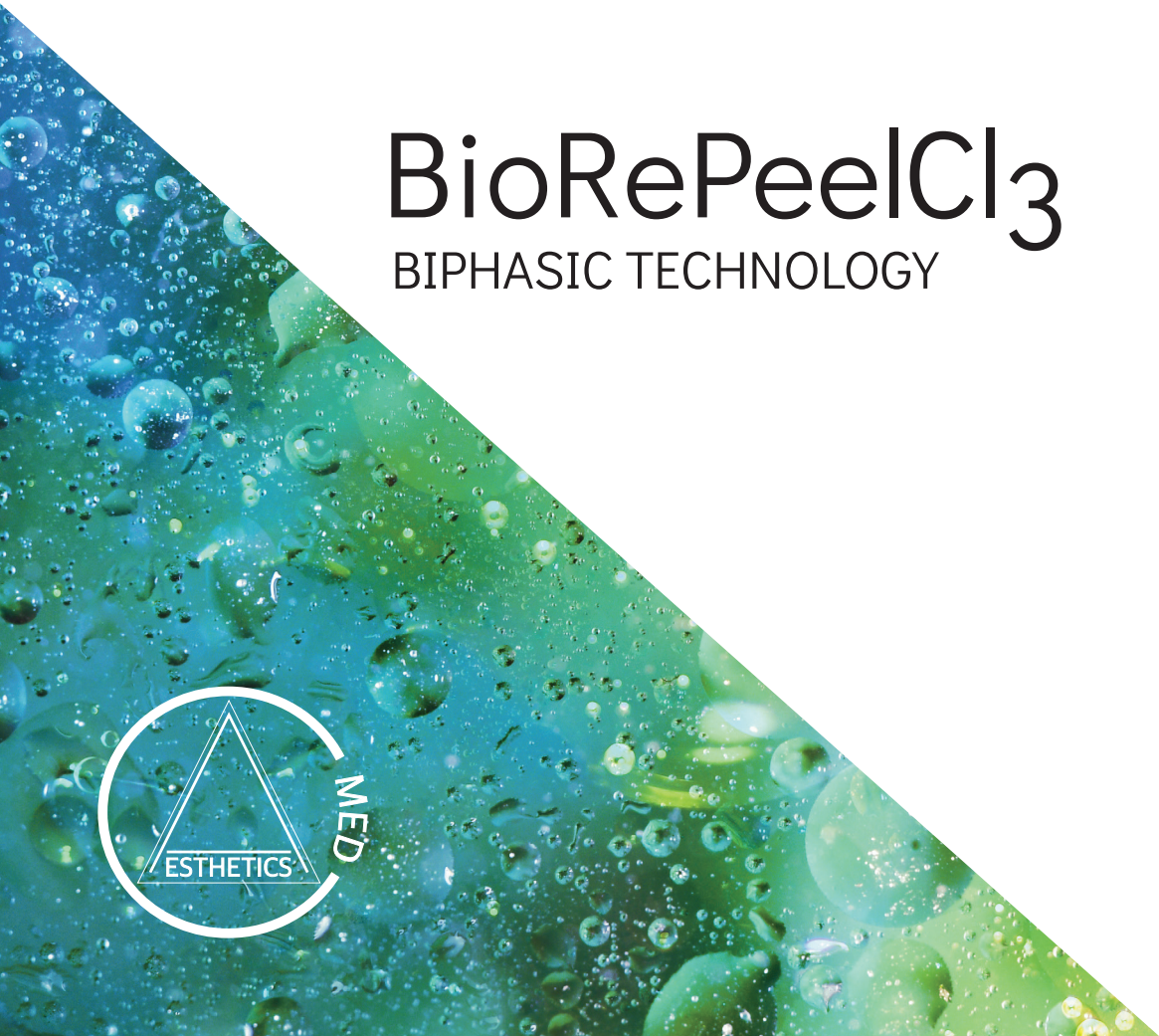

Biostimulation Revitalization Peeling

BioRePeelCl$_3$ is an innovative, biphasic, patented medical device with a biostimulating, revitalizing and peeling action.
Thanks to exfoliating and keratolytic action, it promotes cell turnover by removing the outer corneal layer cells with selective action on corneocytes' junction bonds. The active ingredients' dermofunctionality contributes to a profound stimulation of skin replication and biosynthetic processes, essential to counter the effects of skin aging and oxidative damage induced by UV radiation and pollution.
BioRePeelCl$_3$ is manufactured according to pharmaceutical standards to guarantee the highest quality levels of MADE IN ITALY.

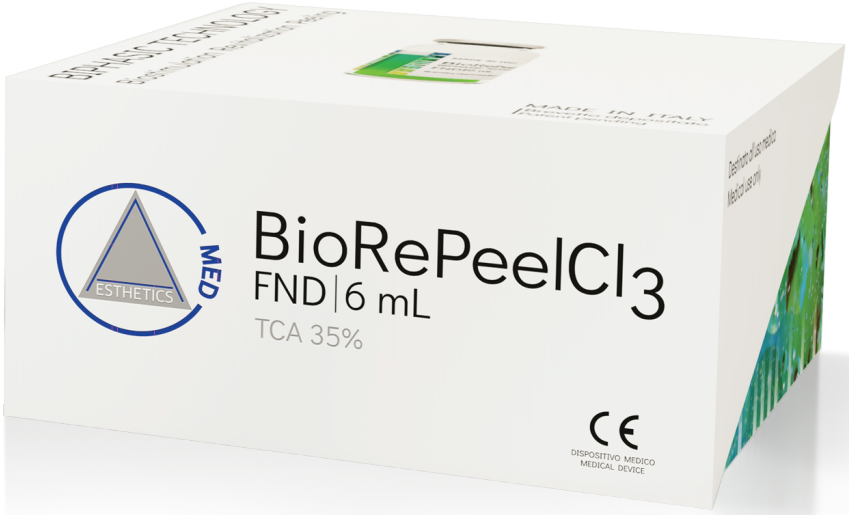



STABILIZING AND PROTECTIVE ACTIONS

It prevents trichloroacetic acid oxidation and stabilizes the formulation separating the aqueous phase from oxidating atmosphere inside the ampoule. It consequently improves performance and effectiveness of the product.

FILMING AND VEHICULATING ACTIONS

It reduces trans-epidermic water evaporation and assists active principles vehiculation.

RIGENERATING ACTION

It restores and maintains cutaneous hydro-lipidic film.

MOISTURIZING ACTION

It increases hydration through direct and indirect mechanisms.



HYDROPHILIC PHASE

BIOSTIMULATING ACTION

It stimulates fibroblasts' activity and activates cutaneous proliferating and biosynthetic processes.

REVITALIZING ACTION

It increases cellular turnover and production of collagen, elastin and hyaluronic acid.

KERATOLYTIC AND PEELING ACTIONS

It removes outer cells of stratum corneum by breaking corneocytes' adesion bonds. This process promotes active principles penetration.

ANTI-AGING AND ANTIOXIDANT ACTIONS

It reduces fine lines by neutralization of free radicals.

ANTI ACNE ENLARGED PORES COMEDONES

It helps battling acne, it eliminates comedones and it has an astringent effect on enlarged pores.



FND|6 mL

FACE ▪ NECK ▪ DÉCOLLETÉ

TRICHLOROACETIC ACID 35%

Lactobionic Acid
Salicylic Acid
Tartaric Acid Citric Acid
Glycine Proline Hydroxyproline Arginine
Vitamin C Vitamin B2
gamma-Aminobutyrric acid
Squalane
Isopropil miristate

PROTOCOL OF USE: Treat the interested part every 7-12 days for 4-6 sessions depending on medical evaluation and treatment targets. Keep the products 3-5 minutes according to patient's perception.

FORMAT:
5 AMPOULES 6 mL

BODY | 12 mL

HANDS ▪ ELBOWS ▪ GLUTES ▪ KNEES ▪ FEET

TRICHLOROACETIC ACID 50%

Lactobionic Acid
Salicylic Acid
Tartaric Acid Citric Acid
Glycine Proline Hydroxyproline Arginine
Vitamin C Vitamin B2
gamma-Aminobutyrric acid
Squalane
Isopropil miristate

PROTOCOL OF USE: Treat the interested part every 10-15 days for 4-6 sessions depending on medical evaluation and treatment targets. Keep the products 3-10 minutes according to patient's perception.

FORMAT:
3 AMPOULES 12 mL



APPLICATION



BioRePeelCl_3 is a biphasic sytem, blue phase is lipophilic and yellow phase is hydrophilic.



Shake the ampoule for a few seconds until colour turns green for a temporary emulsion. Using a syringe, prelate the quantity of product for the treatment.



Apply the product with nitrile gloves until it is completely absorbed. After the laying down-time, clean the zone with a gauze soaked with water.

Emulsion separates in a few minutes and biphasic system is restored. It is necessary to shake the ampoule before every use.



ADVANTAGES

NO NEEDLE	NO OXIDATION DARKENING
NO FROST	COLD STORAGE NOT NECESSARY
INDOLOR	FOR THE WHOLE BODY
FAST	FOR ALL AGES
MULTIFUNCTIONAL	FOR EVERY SEASON
IMMEDIATE RESULT	NO PHOTOSENSIBILIZATION
NO POST-TREATMENT SCARS	NO AMMONIA



ASSOCIATIONS

You can use BioRePeelCl$_3$ also in synergy with other treatments to improve their performance and result:

TREATMENT	WHEN
INJECTABLE TREATMENTS	BEFORE
PHOTOSTIMULATION	BEFORE
RADIOFREQUENCY	AFTER
TOPICAL TREATMENT	BEFORE

MED ESTHETICS

MADE IN ITALY

Patent pending

FEEL FREE TO CONTACT US:
commerciale@cmedaesthetics.com

VISIT OUR WEBSITE:
www.cmedaesthetics.com

CE

MEDICAL DEVICE

CMed Aesthetics S.r.l.
Via Cesare Studiati 4, Pisa, Italy
www.cmedaesthetics.com
commerciale@cmedaesthetics.com

